




**Cordless Neck & Back
Shiatsu Massager with Heat (AM103)**

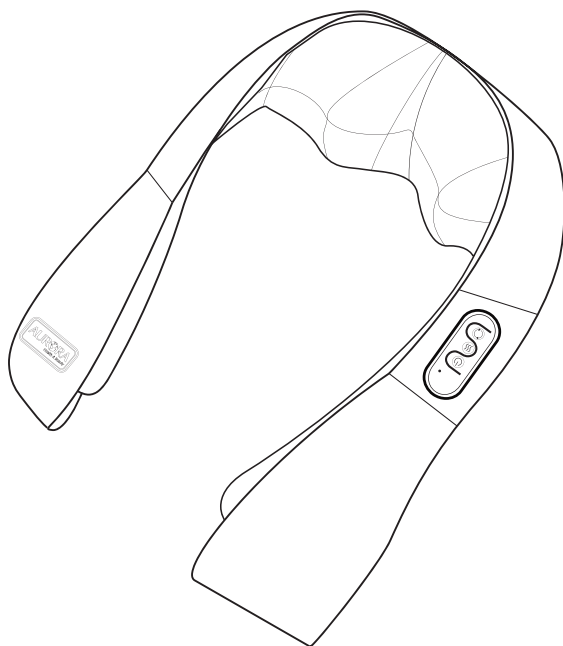

AM103

CAUTION

- TO REDUCE THE RISK OF ELECTRIC SHOCK, DO NOT REMOVE COVER. THERE ARE NO SERVICEABLE PARTS INSIDE.
- TO REDUCE THE RISK OF FIRE OR ELECTRIC SHOCK, DO NOT EXPOSE THIS UNIT TO RAIN OR MOISTURE.

IMPORTANT

SAFETY INSTRUCTIONS

All of the safety and operating instructions should be read, adhered to and followed before the unit is operated.

DANGER

To reduce the risk of electric shock, burns, fire or injury:

1. Always unplug the product from the electrical outlet immediately after using and before cleaning.
2. Do not use while bathing or in a shower.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for product that has fallen into water. Unplug it immediately.

WARNING

1. Do not leave product unattended when plugged in. Unplug from outlet when not in use.
2. This product should not be used on or by children.
3. Do not use attachments not recommended by manufacturer.

4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
5. Keep cord away from heated surfaces.
6. Never block the air openings of the product or place it on a surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint and hair.
7. Never use while sleeping.
8. Do not use outdoors.
9. Never drop or insert objects into any openings.
10. Do not operate where aerosol products are being used or where oxygen is being administered.
11. Do not use on or near eyes or other highly sensitive areas.
12. Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury.
13. Do not carry this product by its power cord or use the cord as a handle.
14. Do not use directly on one concentrated area of the body. It is meant to be used in a continuous motion, constantly changing the massage area.
15. To disconnect, turn all controls to the off position, and then remove the plug from the outlet.
16. Use this appliance only for its intended use as described in this manual. No serviceable parts inside.
17. This product is designed for external body use only.
18. For adult use only.

CAUTION

- Consult your doctor before using this product.
- Individuals with pacemakers should not use this product.
- Never leave the product operating or unattended, especially if children are present.
- Never use the product on open skin eruptions or on swollen or inflamed areas.
- Do not operate after taking painkillers, sedatives or alcohol. Consult your doctor if under these medications.
- Do not operate for more than 20 minutes continuously. It is recommended to turn off the product and allow it to rest (cool) for 40 minutes before re-using. This will prolong the life of the product.
- The product is designed as a non-professional massager to soothe tired and aching muscles. It is not intended to be a substitute for medical attention. If your ailments and symptoms persist or worsen after use, contact your doctor immediately.

LI-ION BATTERY WARNINGS

This device contains a Lithium-ion battery. The battery must be recycled or disposed properly. Contact your local waste management office for information on battery recycling or disposal.

SPECIFICATIONS

Adapter: Input:100-240V~50/60Hz
Output:14V 24W
Battery Type: Li-ion 2500mah
11.1V

Battery Life: Approx. 2.5 hours
(without heat)



Battery Charge Time:
Approx. 4 hours

CARE AND MAINTENANCE

1. Unplug the massager from the electrical outlet before cleaning.
2. Outer cover spot clean only.
3. Store in a cool, dry place.

PARTS LIST

1. DC In Jack (rear)
2. Forward/Reverse Button
3. Heat Button
4. Power Button /High speed/Low speed
5. LED Indicator for charging
6. AC Adapter

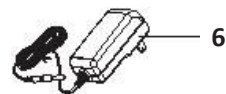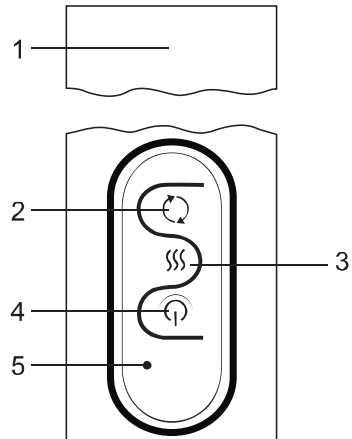

CHARGING

Fully charge the massager before first use.

1. Plug the AC ADAPTER into the DC IN JACK. Plug the AC ADAPTER into an electrical outlet.
2. The CHARGING LED INDICATOR will flash slowly to indicate charging.
3. The CHARGING LED INDICATOR will alight all through when the massager is fully charged.
4. Approximate charging time is 240 minutes (4 Hours).
5. The CHARGING LED INDICATOR will flash red when the battery needs to be charged.

OPERATION

1. Press the POWER BUTTON to turn the massager on with high speed.
Press the POWER BUTTON again, to change the massage to low speed.
Press the POWER BUTTON thirdly, to turn the massage off.
2. Press the FORWARD/REVERSE BUTTON to change the massaging nodes direction.
3. Press the HEAT BUTTON to turn heat on. Press the HEAT BUTTON again to turn heat off.

Note: Massager automatically shuts off after 20 minutes of continuous use.

Warranty & Customer Support

For questions, assistance, support, or warranty claims, please contact us at email address below along with information on where you purchased the product. A copy of your sales receipt/invoice should always be included for proof of purchase. Send all inquiries to:
service@aurorahealthbeauty.com

This product comes with a 1 Year Limited Warranty for more details on our warranty please check the Aurora Health and Beauty website:
www.aurorahealthbeauty.com

Thank you for purchasing an Aurora Health and Beauty product.